



Sycle Noah Sync

Installation and Upgrade Guide



April 17, 2023
Revision 1.0

Copyright and Trademark

Copyright © 2023 Sycle, LLC. All rights reserved. Sycle® is a registered trademark of Sycle, LLC. Other names may be trademarks of their respective holders. Reference to these names is not intended to imply sponsorship, affiliation, or endorsement of Sycle, LLC.

Noah® is a registered trademark of HIMSA II K/S (HIMSA) in the USA. Outside the USA, Noah™ is a trademark of HIMSA II K/S. The Noah software and all intellectual rights pertaining to the Noah software are owned by HIMSA II K/S. Use of the Noah software is subject to the Noah End User License Agreement (EULA) and the terms of use contained in the EULA.

Disclaimer

The information provided by Sycle, LLC ("we," "us," or "our") in this document (the "Document") is subject to change without notice. We may make improvements and/or changes at any time to the product(s) and/or program(s) described in the Document. While all information in the Document is provided in good faith, we provide this information with no representation or warranty of any kind, express or implied.

The Document can not and does not contain medical/health advice. Any examples of medical/health information within the Document are provided with the intention of demonstrating the function of a product or service and are not intended to be seen as accurate or representative medical information.

The Document may contain links to materials or websites containing content belonging to or originating from third parties. We do not warrant, endorse, guarantee, or assume responsibility for the accuracy or reliability of any information or service offered by third-party material or websites linked within the Document.

The Document and the information provided in the Document are not authorized by, endorsed by, sponsored by, affiliated with, or otherwise approved by HIMSA II K/S.

Table of Contents

1. Overview	1
1.1. Installation Summary	1
1.2. Upgrade Summary	1
2. Obtaining an Upload Link and a Conversion License	2
3. Installing Sycle Noah Sync	3
3.1. Installation Requirements	3
3.2. Migrating Noah System Data	4
3.2.1. Creating a Noah System Data Backup	4
3.2.2. Exporting and Uploading Noah System Data	7
3.2.3. Converting the Noah License	13
3.3. Running the Installation	15
4. Upgrading Sycle Noah Sync	19
4.1. Sycle Noah Sync (SNS) Upgrade Requirements	19
4.2. Backing Up the Noah Directory	20
4.3. Running the Upgrade	22
A. Glossary	27
B. Document Revision History	28

1. Overview

This guide covers the preparation, data migration, and installation instructions for moving your practice to [Sycle Noah Sync \(SNS\)](#).

This process and guide are structured accordingly for installing or upgrading SNS. Summaries of the instructions are below.

1.1. Installation Summary

This is a summary of the workflow and instructions provided for a first-time installation of SNS. This summary (and the full chapter) assumes that you have an existing Noah installation on either a standalone Noah database (on a single PC) or using a shared office database.

In the [Obtaining an Upload Link and a Conversion License](#) chapter:

- Days (preferably about a week) prior to migrating your Noah data and installing Sycle Noah Sync, you contact Sycle's Customer Experience Team (at **888.881.7925** or <support@sycle.net>) to obtain
 - an upload link for your exported Noah patient data, and
 - a Noah Conversion key (used to convert your Noah System license to a Noah Engine license)

In the [Installing Sycle Noah Sync](#) chapter, you are instructed to:

- Review the requirements before proceeding with the installation ([Installation Requirements, page 3](#))
- Back up your Noah System data ([Creating a Noah System Data Backup, page 4](#))
- Export and upload your Noah System data file to Sycle ([Exporting and Uploading Noah System Data, page 7](#))
- Convert your Noah System license to a Noah Engine license ([Converting the Noah License, page 13](#))
- Install the software ([Running the Installation, page 15](#))

1.2. Upgrade Summary

This is a summary of the workflow and instructions provided for upgrading an existing SNS installation.

In the [Upgrading Sycle Noah Sync](#) chapter, you are instructed to:

- Review the requirements before proceeding with the upgrade ([Sycle Noah Sync \(SNS\) Upgrade Requirements, page 19](#))
- Back up your Noah directories ([Backing Up the Noah Directory, page 20](#))
- Install the upgrade ([Running the Upgrade, page 22](#))

2. Obtaining an Upload Link and a Conversion License

Prior to migrating your Noah data and installing Sycle Noah Sync, you need to contact the [Sycle Customer Experience Team](#). Ideally, this should be done some days (preferably a week) before you perform the migration.

CE will provide you with a link that you will use to upload your exported Noah System data (in the [Exporting and Uploading Noah System Data](#) chapter). They will also provide a Noah Conversion key that you will use to convert your Noah System license to a Noah Engine license (in the [Converting the Noah License](#) chapter).

You can contact the Customer Experience team by calling 888.881.7925 or by emailing <support@sycle.net>.



If you already have SNS installed, do *not* proceed with the instructions in this chapter. Instead, refer to [Upgrading Sycle Noah Sync, page 19](#).

If you choose to call the CE team at 888.881.7925:

- ☐ Be ready to provide CE with an email address that you can access on the PC containing your Noah installation.
- ☐ Have the registration number (a 16 digit alpha-numeric code) associated with your current Noah license ready. Your registration code can be found on the Proof-of-Purchase certificate.
- ☐ Inform CE that you are performing a migration from Noah (standalone) to Sycle Noah sync.
- ☐ Provide CE with a timeline for the migration.

If you email the CE team at <support@sycle.net>, do the following:

- ☐ Ensure you can access the email address used to send the message to CE on the PC containing your Noah installation.
- ☐ Indicate that you are performing a migration from Noah (standalone) to Sycle Noah sync in the title/body of the email.
- ☐ Provide a timeline for the installation and migration.
- ☐ Provide the registration number (a 16 digit alpha-numeric code) associated with your current Noah license in the email message. Your registration code can be found on the Proof-of-Purchase certificate.
- ☐ Confirm the email address used to send the message is the correct email address to send your conversion key and upload link.

When you have the upload link and conversion key, you are ready to perform the migration and installation. Proceed to [Installing Sycle Noah Sync, page 3](#).

3. Installing Cycle Noah Sync

Installing Cycle Noah Sync (SNS) is a multi-step process that involves backing up your Noah System data, exporting and uploading your Noah System data file, converting your Noah System license to a Noah Engine license, and installing the software.



If you already have SNS installed, do *not* proceed with the instructions in this chapter. Instead, refer to [Upgrading Cycle Noah Sync, page 19](#).

3.1. Installation Requirements

Before you proceed with installing SNS, ensure that the following requirements are met. You may require assistance from your system administrator or IT department.

- ☐ You will not need to make any changes to NOAH patient data in the next 24–36 hours. Cycle recommends that you perform this migration after business hours and/or on a Friday to minimize any operation impact.
- ☐ If you have any third-party Noah modules installed, you have consulted your module manufacturers and confirmed that those modules are compatible with the version of Noah incorporated in SNS.
- ☐ Your PC is running one of the operating systems supported by SNS and Noah System 4.
 - Windows 11 – Excluding Windows 11 on Arm and Insider Preview builds.
 - Windows 10 – Excluding Long Term Service Channel/Branch (LTSC/LTSB) releases, Windows 10 Mobile OS, Windows 10 in S mode, Windows 10 on Arm, and Insider Preview builds.
 - Windows Server 2022/2019/2016 – Excluding Windows Server Core and Nano Server.
- ☐ Your PC meets the other minimum system requirements for Noah System 4¹.
- ☐ Your Windows operating system has the latest Windows Updates and Service Packs applied.
- ☐ You know the administrator username and password for your PC. The installer must be run as an administrator account².
- ☐ (If applicable) Your system administrator or IT department are aware of the installation timeline and are available to offer assistance.
- ☐ You have obtained and have access to the upload link provided by the [Cycle Customer Experience Team](#) (see [Obtaining an Upload Link and a Conversion License, page 2](#)).
- ☐ You have obtained and have access to the Noah Conversion key (provided by the [Cycle Customer Experience Team](#)) that you will use to convert your Noah System license to a Noah Engine license (see [Obtaining an Upload Link and a Conversion License, page 2](#)).

¹These requirements are listed in the [Minimum system requirements for Noah System 4](#) article on <http://www.himsa.com>.

²The Administrator account must be set up with a valid password. Running the installer with an Administrator account without a password will fail.

3.2. Migrating Noah System Data

The Noah System's patient data must be backed up (for redundancy), exported, uploaded to Cycle, and your Noah System license converted to a Noah Engine license before you install Cycle Noah Sync.

3.2.1. Creating a Noah System Data Backup

The patient data on your Noah System must be backed up. This ensures you can restore your Noah database if there are problems with the migration.

These instructions assume that the default database backup location is configured. You can verify your configured backup location by in the Noah Console by navigating to *Database Tools* → *Backup Scheduler* and reviewing the *Backup Location* → **Location** field.

If you are using a custom backup location outside the `C:\ProgramData\HIMSA*` directory (such as an external storage device or a network share), follow only [Step 1](#) and [Step 2](#).



These steps must be performed immediately before installing Cycle Noah Sync, and at a time when no changes will be made to patient data for the next 24–36 hours (see [Installation Requirements, page 3](#)).



This guide assumes that you are migrating either a standalone Noah database (on a single PC) or using a shared office database. If using a shared office database, ensure all steps are performed on the PC running the shared office database (i.e. your Noah Server).

1. Open the **Noah Console**.




The **Noah Console** can be opened from the Noah software (*File* → **Noah Console**) or from Windows (*Programs* → *HIMSA* → **Noah Console**)

2. Navigate to and click *Noah Console Tools* → *Database Administration* → *Tools* → **Backup**

3. In Windows Explorer, navigate to `C:\ProgramData\HIMSA`.



If you can not see the `ProgramData` directory in `C:\`, enter the full path directly in to your Windows Explorer address bar and press **Enter**.

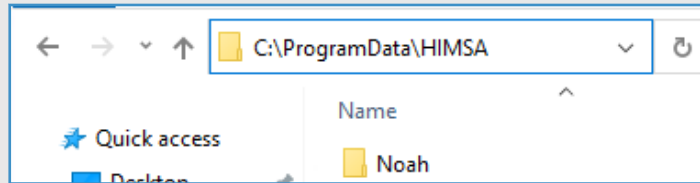


If the above directory does not exist, navigate instead to `C:\Program Files (x86)\HIMSA`.

4. Copy the Noah directory from C:\ProgramData\HIMSA to a secondary location. This example uses the local PC's Desktop.




If the above directory does not exist, copy the Noah 4 directory from C:\Program Files (x86)\HIMSA to the secondary location using the same steps.


- a. Click the Noah directory and press **CTRL+C** to copy the directory.

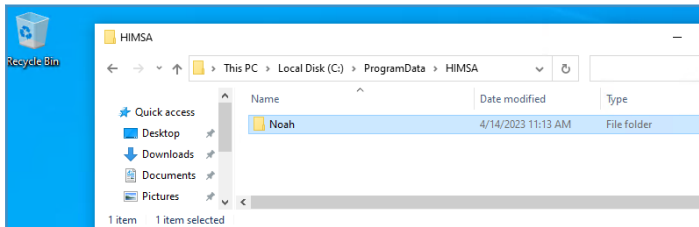

- b. Navigate to your secondary location. In this example, we navigate to the Desktop (⌘ Win+D).

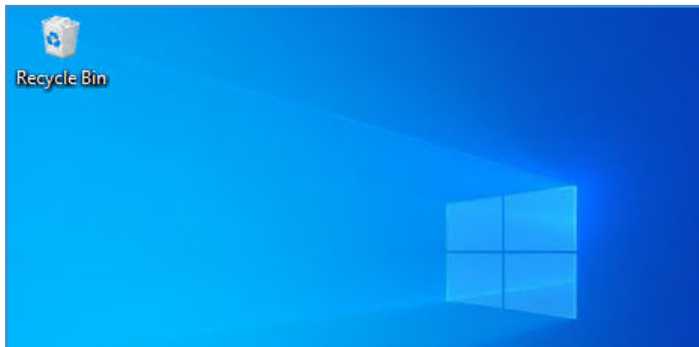

- c. Press **CTRL+V** to paste the directory to your secondary location.

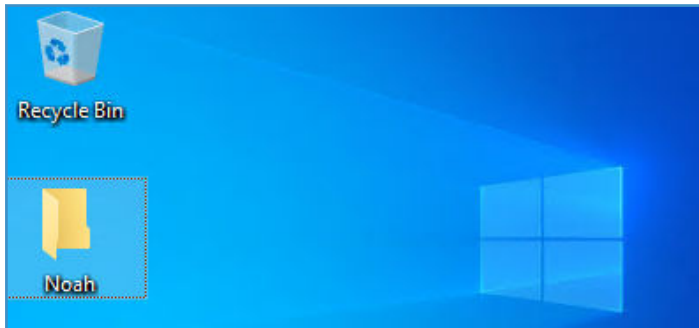

3.2.2. Exporting and Uploading Noah System Data

After you have backed up your Noah System data ([Creating a Noah System Data Backup, page 4](#)), you must export your patient data and upload it to Sycle.




These steps must be performed immediately before installing Sycle Noah Sync, and at a time when no changes will be made to patient data for the next 24–36 hours (see [Installation Requirements, page 3](#)).




This guide assumes that you are migrating either a standalone Noah database (on a single PC) or using a shared office database. If using a shared office database, ensure all steps are performed on the PC running the shared office database (i.e. your Noah Server).

Exporting Noah System Data

1. In your Noah software, navigate to and click *File* → **Export Patients**.

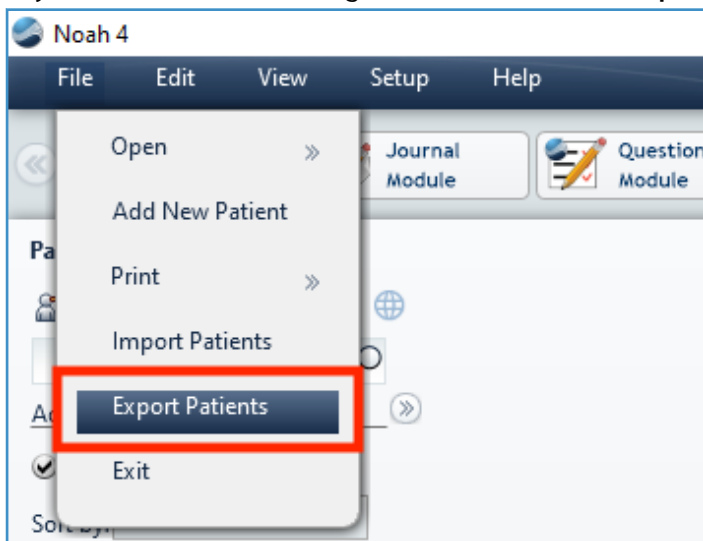

2. In the *Export Options* window, do the following:

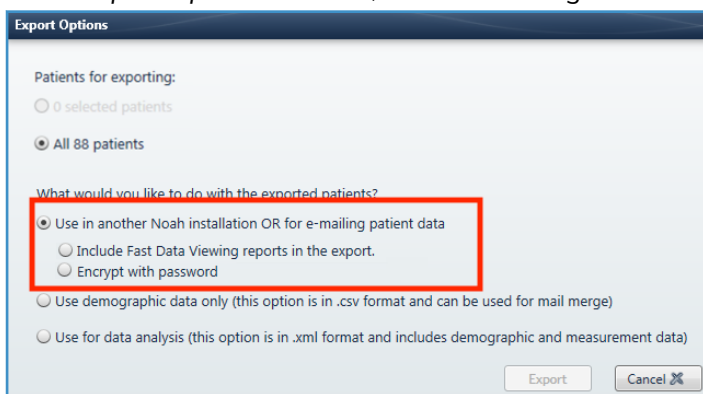

- a. Select *Under Patients for exporting* → **All patients**.
- b. Select *What would you like to do with exported patients?* → **Use in another Noah installation OR for e-mailing patient data**.
- c. Under **Use in another Noah installation OR for e-mailing patient data**, deselect the **Include Fast Data Viewing reports in the export** and **Encrypt with password** options.




These options are enabled by default. Ensure these fields are *not* enabled before proceeding. Disabling these options is required for Noah to export the patient data in the .nhax format required for Sycle to import your data to Sycle Noah Sync.



These fields will not appear on versions of Noah lower than Noah 4.9.

- d. Click **OK** or **Export**.

3. In the *Export Note* window, enter the *Intended Recipient* as **Sycle, LLC** and the *Purpose of Export* as **Migration**. Click **Save**.

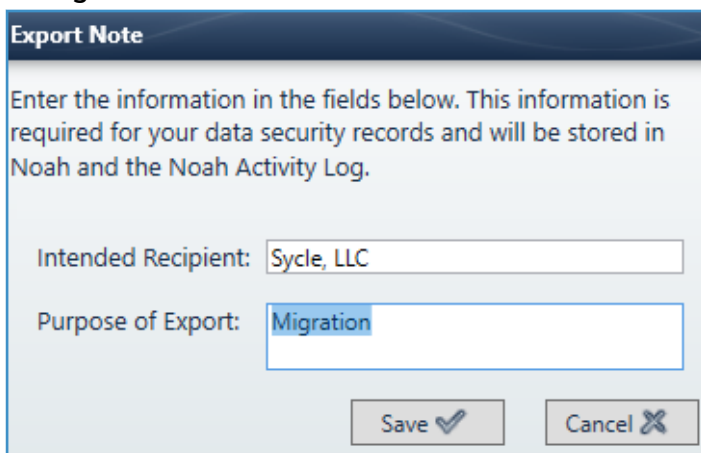

Export Note

Enter the information in the fields below. This information is required for your data security records and will be stored in Noah and the Noah Activity Log.

Intended Recipient:

Purpose of Export:




This window will not appear on versions of Noah lower than Noah 4.10.

4. In the *Save As* window, navigate to a location where you will save the file, then do the following:

a. In the *File Name* field, enter the name of the file in the following format:

PCXXXX-Clinic Name-Date

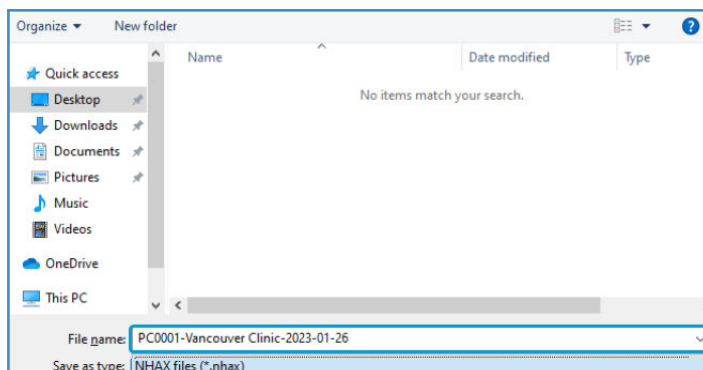


Dates should be in ISO 8601 format, e.g. YYYY-MM-DD.



An example file name would be PC0001-Vancouver Clinic-2023-01-26

b. In the *Save as Type* field, select **NHAX files (*.nhax)**.

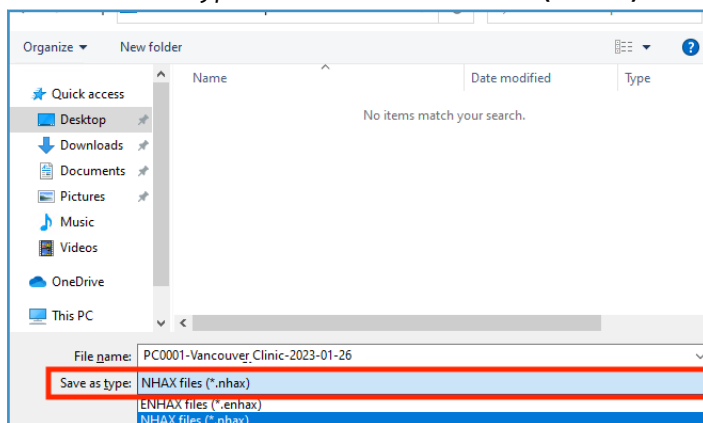

c. Click **Save**.

Uploading Noah Data to Sharefile

1. Open the email provided by the Customer Experience team.




This email was sent as a part of the preparation steps performed in the [Obtaining an Upload Link and a Conversion License](#) chapter.

2. Click the link in the email from the Customer Experience team.




Never click any links in *unsolicited* emails. If in any doubt about an email message that seems to originate from Sycle, contact the *.Sycle Customer Experience Team*.


3. Verify the domain (`securesycle.sharefile.com`) in your browser's address bar and that the CE team member's name is shown on the page.

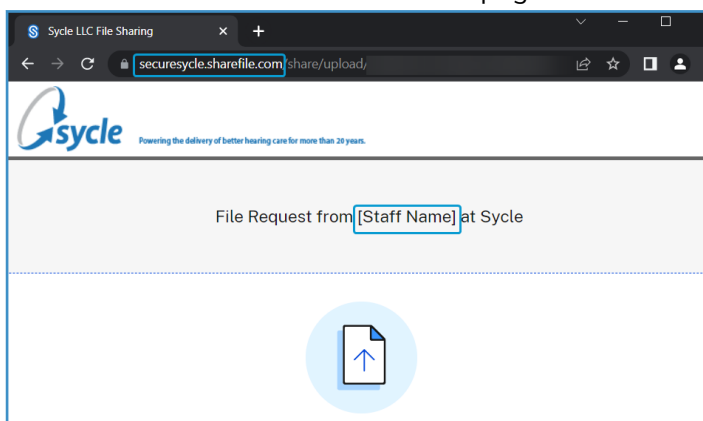


4. Click **Browse Files**.

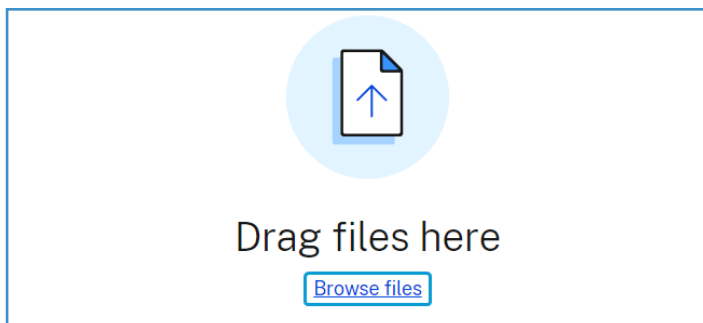


5. In the Open window, navigate to and select the PCXXXX-Clinic Name-Date file you exported in [Step 4](#) of the [Exporting Noah System Data](#) instructions chapter and click **Open**.

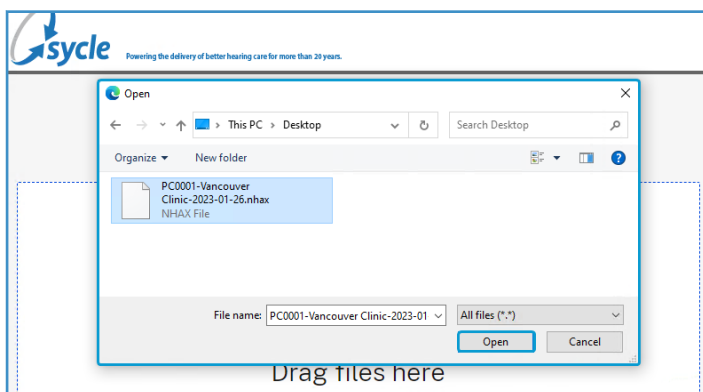

6. Verify the information on screen, then click **Upload**.

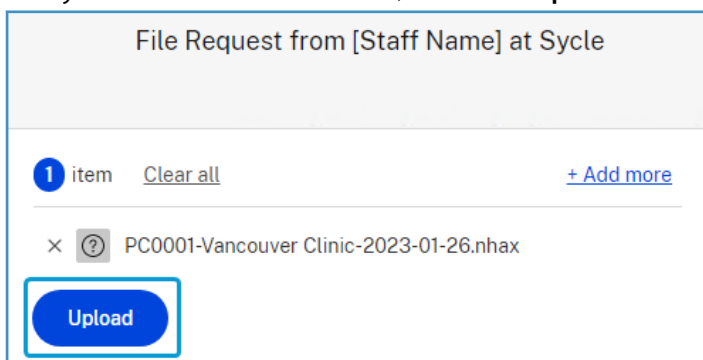

7. Wait for the file upload to complete.

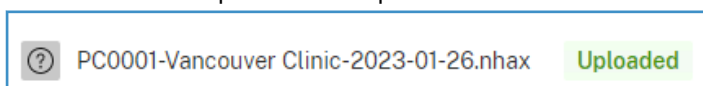

After the file has been uploaded, contact the [Sycle Customer Experience Team](#) team by calling 888.881.7925 or by emailing support@sycle.net and verify that the file upload was successful. If the file upload was successful, proceed to the [Converting the Noah License, page 13](#) chapter.


3.2.3. Converting the Noah License

You will now use your Noah Conversion key to convert your Noah System license to a Noah Engine license.

You will need both the Noah Conversion key provided by the [Sycle Customer Experience Team](#) (see [Obtaining an Upload Link and a Conversion License, page 2](#)) and the registration number associated with your current Noah license³.


Converting to and Downloading the Noah Engine License

1. Navigate to the *Noah License Registration Page* at registration.himsa.dk




You can also open the *Noah License Registration Page* by navigating to himsa.com, then selecting *Quick Links* → **Register your Noah license**.


2. In the **License Registration Number** field, enter the **Noah Conversion key** provided by the Customer Experience Team and click **Submit**.

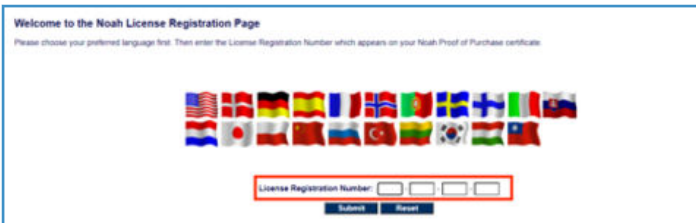

3. Review the *License Agreement* and click **I Agree**.
4. On the *Noah upgrade* page, enter the registration number associated with your current Noah license³ in the **Registration Number** field and click **Submit**.



5. On the *Register your Noah License* page, verify the accuracy of the information in all fields and click **Submit**.



You must contact HIMSA to resolve any inaccuracies on this page.

6. Click **Download** to download the `License.exe` file.

Proceed to [Running the License Conversion File, page 14](#).


³This registration number is the same 16 digit alpha-numeric code you provided to CE in the [Obtaining an Upload Link and a Conversion License, page 2](#) chapter.

Running the License Conversion File


1. Navigate to your Downloads directory.



On most browsers, the easiest way to view your recent downloads is by opening the browser and pressing **CTRL+J** on your keyboard. When the downloads tab or window appears, you can click the appropriate icon or text to open your downloads folder (shown below).

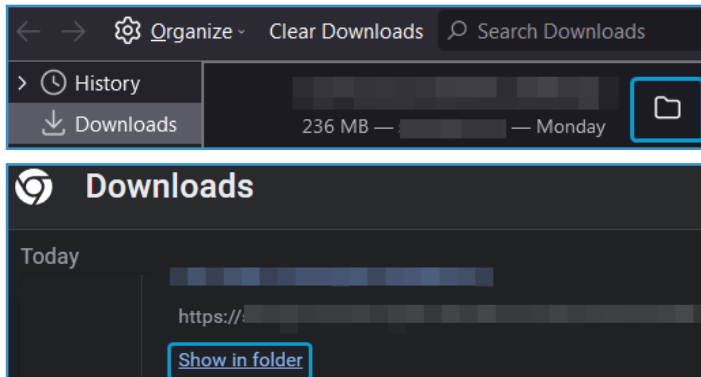

2. Right-click the `License.exe` file and select **Run as Administrator**.

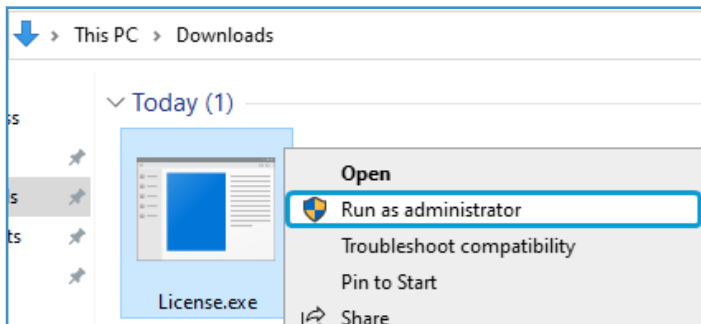

If file name extensions are not visible on your computer, this file will display only as `License` and will not show the `.exe` extension. This does not affect the installation process.

3. If the *User Account Control* window appears, click **Yes**.



This window may not appear on your device. This does not affect the installation process and you can ignore this step.


4. When the installer opens, follow the installation instructions on screen to proceed through the installation.
5. When the installation of the license file is complete, ensure the *Result* window shows "The NOAH license file has been successfully installed" and click **OK**.

3.3. Running the Installation


1. (If applicable) Pause or disable any non-Microsoft antivirus programs or firewalls installed on your PC.
2. Open your web browser and navigate to sycle.com/noah-synch-download.
3. Click the download link on that page to download the SNS installer package.
4. Navigate to your Downloads directory.




On most browsers, the easiest way to view your recent downloads is by opening the browser and pressing **CTRL+J** on your keyboard. When the downloads tab or window appears, you can click the appropriate icon or text to open your downloads folder (shown below).

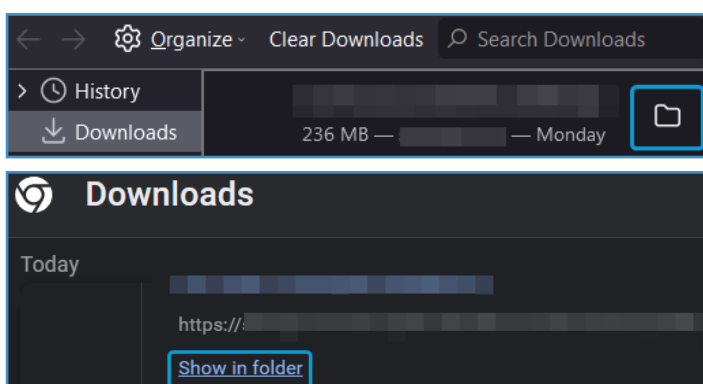

5. In the Downloads directory, **right-click** the downloaded installation package and select **Extract All**.

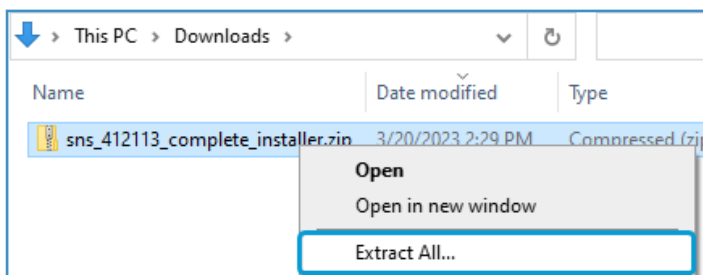


6. Do not change any of the options in the *Extract Compressed (Zipped) Folders* window. Click **Extract**.

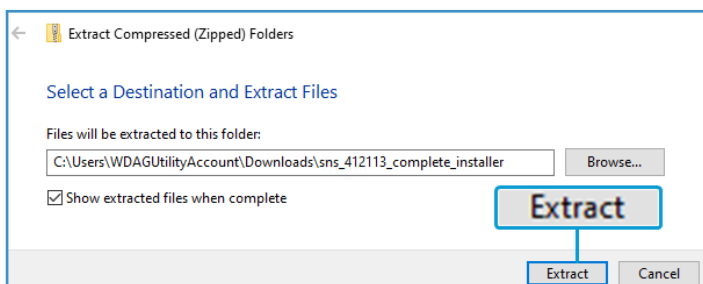


Wait for the files to finish extracting. When the extraction is complete, the location of the extracted files will open in a new Windows Explorer window.

7. **Right-click** the `setup.exe` file and select **Run as Administrator**.

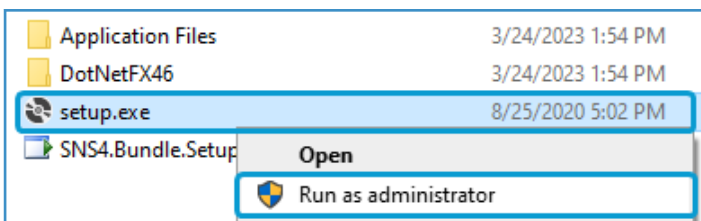


If file name extensions are not visible on your computer, this file will display only as `setup` and will not show the `.exe` extension. This does not affect the installation process.

8. If the *User Account Control* window appears, click **Yes**.

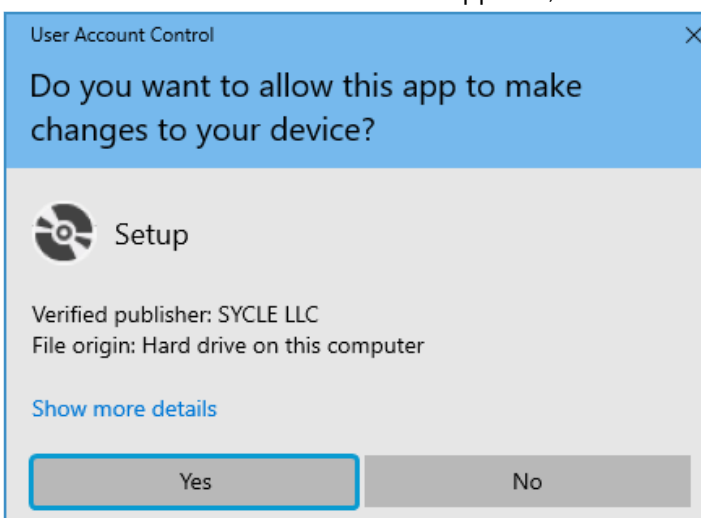


This window may not appear on your device. This does not affect the installation process; simply ignore this step and proceed to [Step 9](#).

9. Wait for the *Launching Application* window to finish verifying the application requirements. Do *not* exit this window manually.

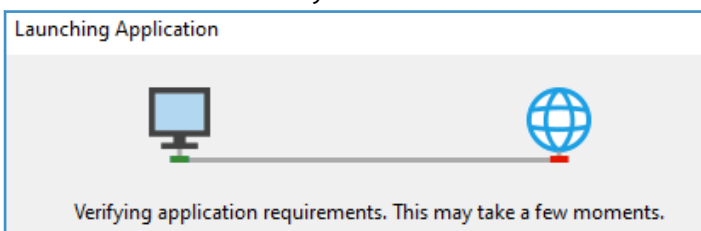


When this step is completed (automatically), the installer closes this window and opens the *Application Run - Security Warning* window.

10. In the *Application Run - Security Warning* window, select **Run**.

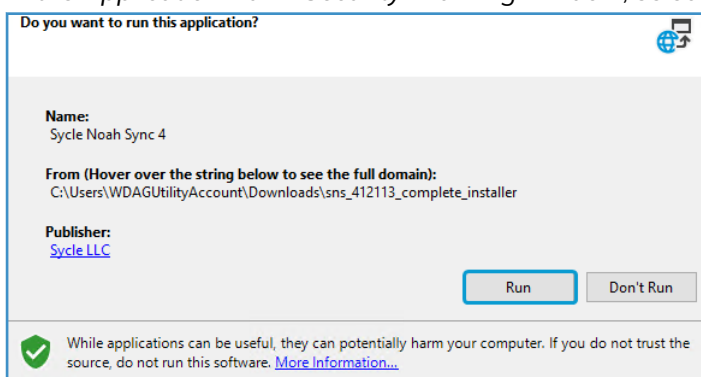

11. Wait for the *Downloading Sycle Noah Sync* window to finish the download process. Do not exit this window manually.

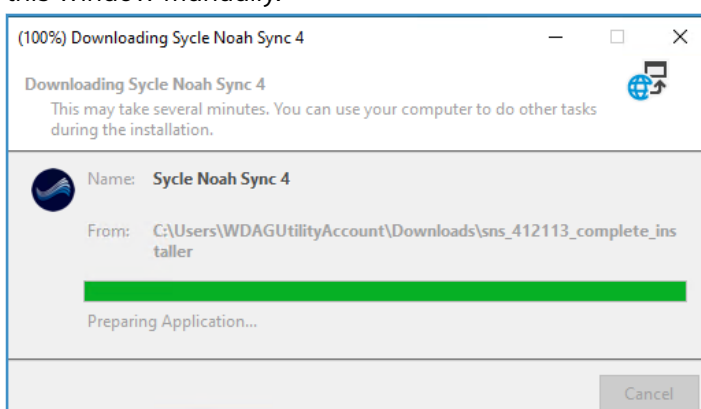


When this step is completed (automatically), the installer closes this window and opens the *Sycle Noah Sync Setup* window.




This window may only be briefly visible before it closes automatically. If you do not see this window, simply proceed to [Step 12](#).

12. Review the contents of the *Sycle Noah Sync Setup* window. If you have fulfilled the conditions and requirements (see [Installation Requirements, page 3](#)), click **Proceed**.

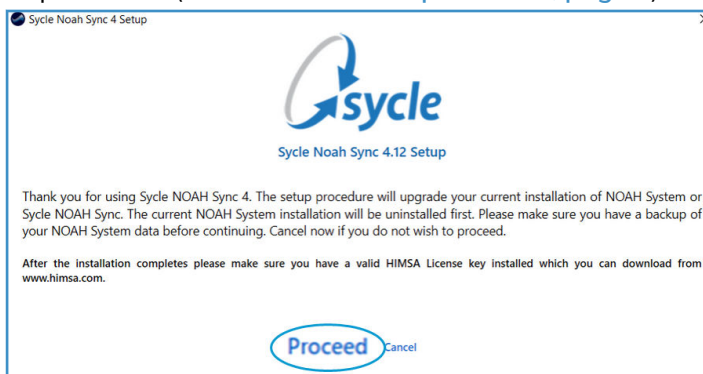

If you have *not* fulfilled the conditions and requirements (see [Installation Requirements, page 3](#)), this is your last chance to cancel the installation process. If you have not reviewed the requirements, click **Cancel** to abort the installation. When all the requirements have been met, you can start the installation again from [Step 7](#).

13. Wait for the installer to complete the installation. This process is automatic. Do not close or exit any of the windows that appear during the installation process.
14. (If applicable) Un-pause or enable any non-Microsoft antivirus programs or firewalls that you disabled in [Step 1](#).
15. When the installer shows **Installation Finished**, click **Close**.

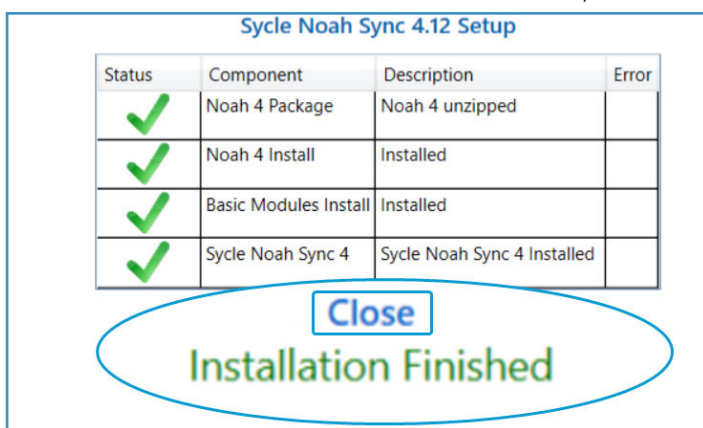

4. Upgrading Cycle Noah Sync

4.1. Cycle Noah Sync (SNS) Upgrade Requirements



These instructions are applicable to SNS 4.9 or greater. If you are upgrading from a version below 4.9, do not run the upgrade. Contact Cycle's Customer Experience Team (at 888.881.7925 or <support@cycle.net>) for assistance.

Ensure that the following requirements are met before you proceed with the upgrade. You may require assistance from your system administrator or IT department.

- ☐ SNS can not be used during the upgrade. The upgrade will take 10–30 minutes depending on your PC and configuration.
- ☐ If you have any third-party Noah modules installed, you have consulted your module manufacturers and confirmed that those modules are compatible with the version of Noah incorporated in the SNS upgrade.
- ☐ You have confirmed that you are upgrading using a tested upgrade path⁴.
- ☐ Your PC is running one of the operating systems supported by SNS and Noah System 4.
 - Windows 11 – Excluding Windows 11 on Arm and Insider Preview builds.
 - Windows 10 – Excluding Long Term Service Channel/Branch (LTSC/LTSB) releases, Windows 10 Mobile OS, Windows 10 in S mode, Windows 10 on Arm, and Insider Preview builds.
 - Windows Server 2022/2019/2016 – Excluding Windows Server Core and Nano Server.
- ☐ Your PC meets the other minimum system requirements for Noah System 4⁵.
- ☐ Your Windows operating system has the latest Windows Updates and Service Packs applied.
- ☐ You know the administrator username and password for your PC. The installer must be run as an administrator account⁶.
- ☐ You, your system administrator, or your IT department are aware of any non-Microsoft antivirus programs or firewalls installed on your PC and are able to disable them for the duration of the installation process.
- ☐ (If applicable) Your system administrator or IT department are aware of the upgrade timeline and are available to offer assistance.

⁴You can review HIMSA's suggested upgrade path in the [Upgrade Path for Unsupported Noah System Versions](#) article, available on www.himsa.com.

⁵These requirements are listed in the [Minimum system requirements for Noah System 4](#) article on <http://www.himsa.com>.

⁶The Administrator account must be set up with a valid password. Running the installer with an Administrator account without a password will fail.


4.2. Backing Up the Noah Directory

You must back up your Noah files prior to running the upgrade. This ensures you can restore your Noah data and modules if there are problems with the upgrade.

1. In Windows Explorer, navigate to `C:\ProgramData\HIMSA`.



If you can not see the `ProgramData` directory in `C:\`, enter the full path directly in to your Windows Explorer address bar and press **Enter**.

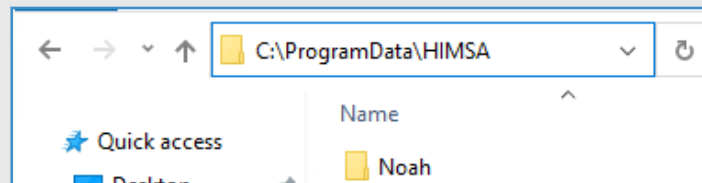


If the above directory does not exist, navigate instead to `C:\Program Files (x86)\HIMSA`.

2. Copy the Noah directory from C:\ProgramData\HIMSA to a secondary location. This example uses the local PC's Desktop.




If the above directory does not exist, copy the Noah 4 directory from C:\Program Files (x86)\HIMSA to the secondary location using the same steps.


- a. Click the Noah directory and press **CTRL+C** to copy the directory.

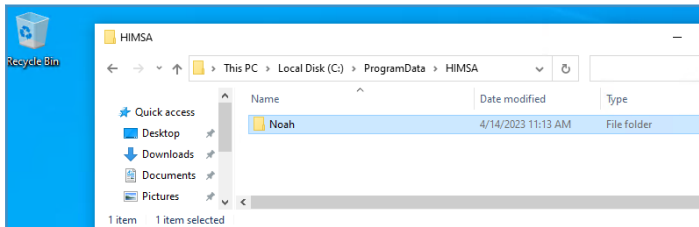

- b. Navigate to your secondary location. In this example, we navigate to the Desktop (⌘ Win+D).

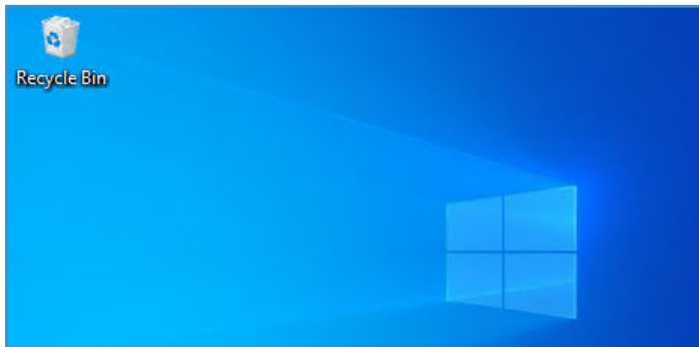

- c. Press **CTRL+V** to paste the directory to your secondary location.

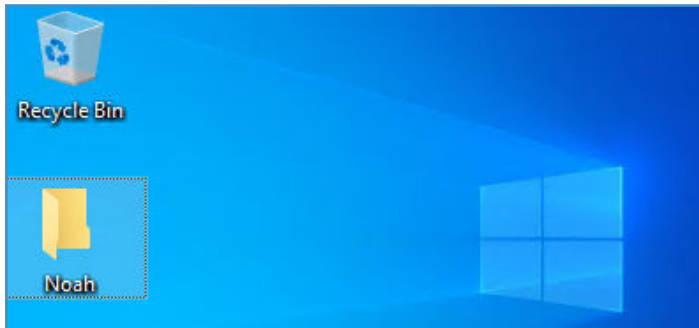

With the backup in place, proceed to [Running the Upgrade, page 22](#).

4.3. Running the Upgrade




These instructions are applicable to SNS 4.9 or greater. If you are upgrading from a version below 4.9, do not run the upgrade. Contact Sycle's Customer Experience Team (at **888.881.7925** or <support@sycle.net>) for assistance.


1. (If applicable) Pause or disable any non-Microsoft antivirus programs or firewalls installed on your PC.
2. Open your web browser and navigate to sycle.com/noah-synch-download.
3. Click the download link on that page to download the SNS upgrade package.
4. Navigate to your Downloads directory.




On most browsers, the easiest way to view your recent downloads is by opening the browser and pressing **CTRL+J** on your keyboard. When the downloads tab or window appears, you can click the appropriate icon or text to open your downloads folder (shown below).

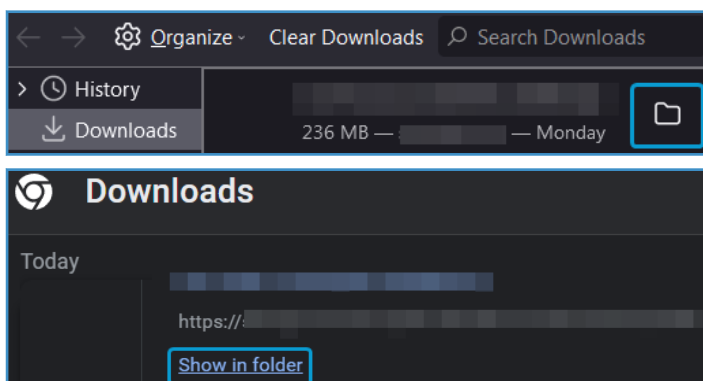

5. In the Downloads directory, **right-click** the downloaded upgrade package and select **Extract All**.

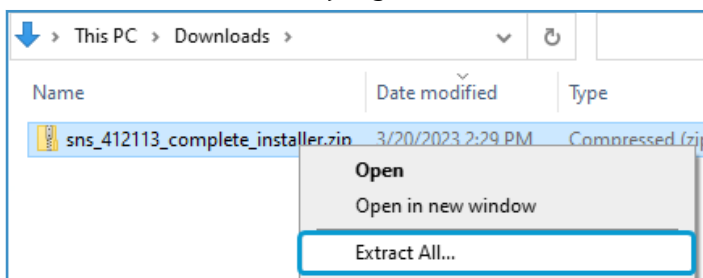


6. Do not change any of the options in the *Extract Compressed (Zipped) Folders* window. Click **Extract**.

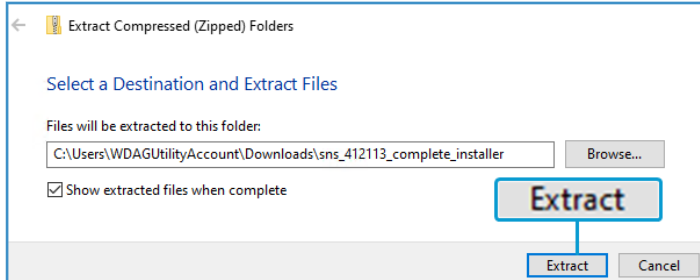


Wait for the files to finish extracting. When the extraction is complete, the location of the extracted files will open in a new Windows Explorer window.

7. **Right-click** the `setup.exe` file and select **Run as Administrator**.

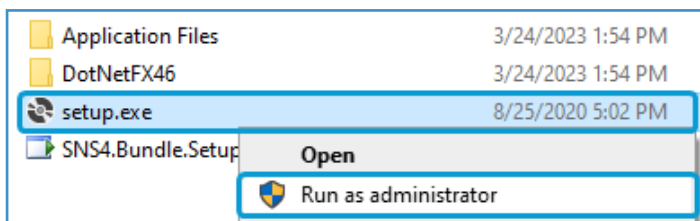


If file name extensions are not visible on your computer, this file will display only as `setup` and will not show the `.exe` extension. This does not affect the upgrade process.

8. If the *User Account Control* window appears, click **Yes**.

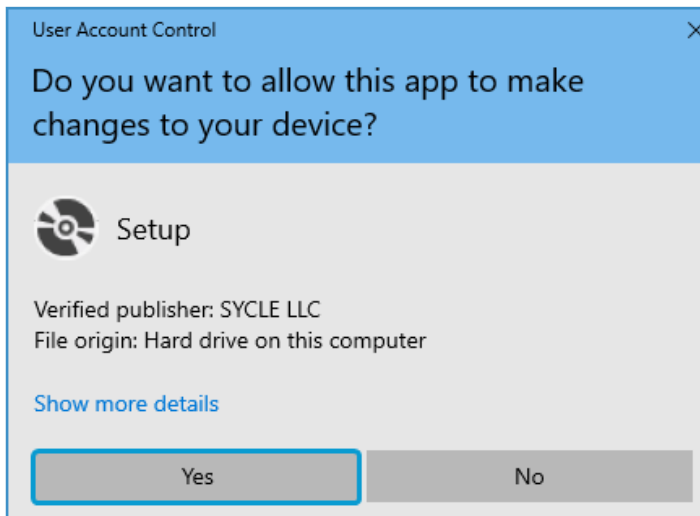

This window may not appear on your device. This does not affect the upgrade process; simply ignore this step and proceed to [Step 9](#).

9. Wait for the *Launching Application* window to finish verifying the application requirements. Do not exit this window manually.

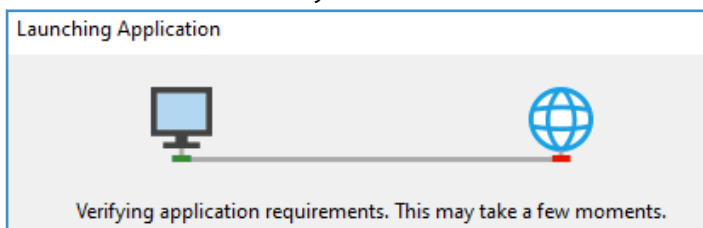


When this step is completed (automatically), the installer closes this window and opens the *Application Run - Security Warning* window.




This window may only be briefly visible before it closes automatically. If you do not see this window, simply proceed to [Step 10](#).

10. In the *Application Run - Security Warning* window, select **Run**.

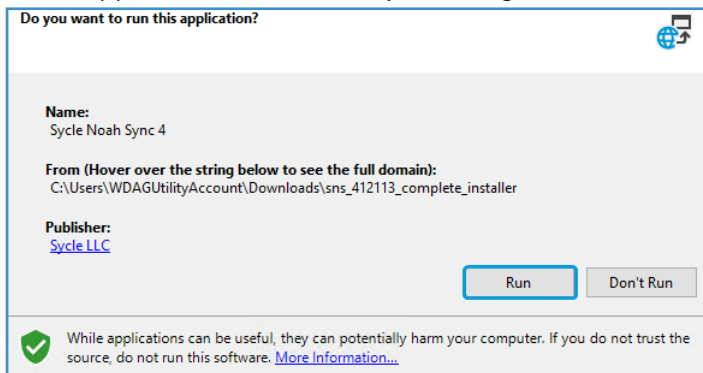

11. Wait for the *Downloading Sycle Noah Sync* window to finish the download process. Do not exit this window manually.

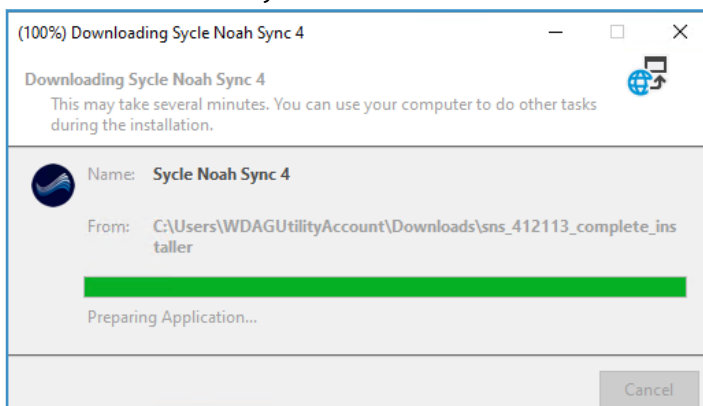


When this step is completed (automatically), the installer closes this window and opens the *Sycle Noah Sync Setup* window.




This window may only be briefly visible before it closes automatically. If you do not see this window, simply proceed to [Step 12](#).

12. Review the contents of the *Sycle Noah Sync Setup* window. If you have fulfilled the conditions and requirements (see [Sycle Noah Sync \(SNS\) Upgrade Requirements, page 19](#)), click **Proceed**.

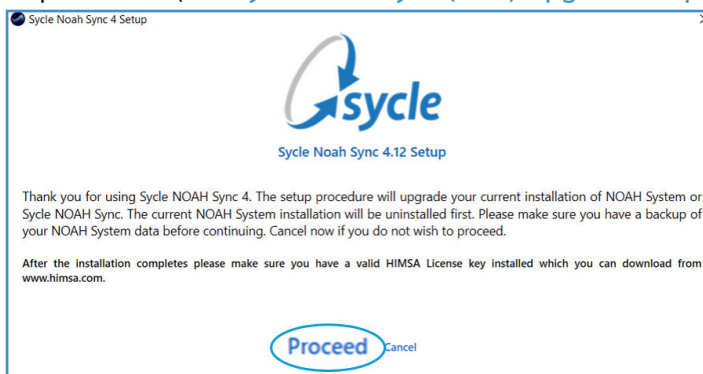


If you have *not* fulfilled the conditions and requirements (see [Sycle Noah Sync \(SNS\) Upgrade Requirements, page 19](#)), this is your last chance to cancel the upgrade process. If you have not reviewed the requirements, click **Cancel** to abort the upgrade. When all the requirements have been met, you can start the upgrade again from [Step 7](#).

13. Wait for the installer to complete the upgrade. This process is automatic.
If you are prompted to close any applications before continuing the upgrade, select **Automatically close applications and attempt to restart them after setup is complete** and click **OK**.

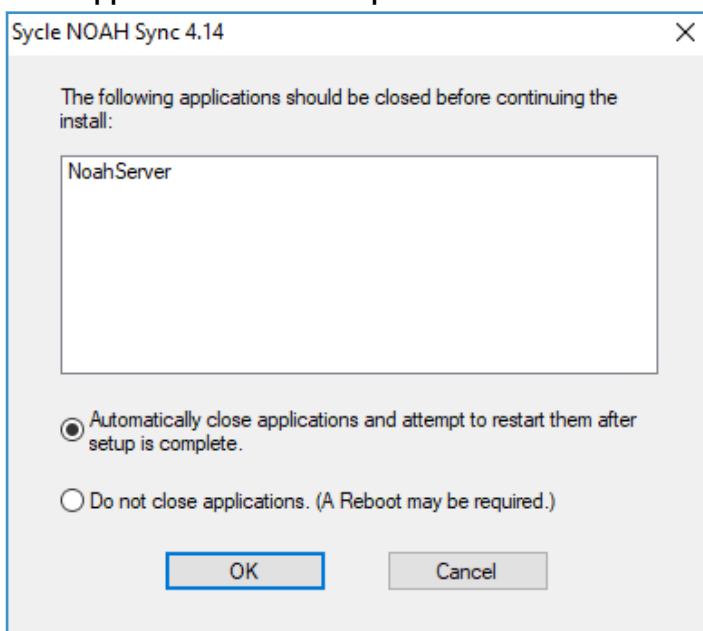

Do *not* close or exit any other windows that appear during the upgrade process.

14. When the installer shows **Installation Finished**, click **Close**.

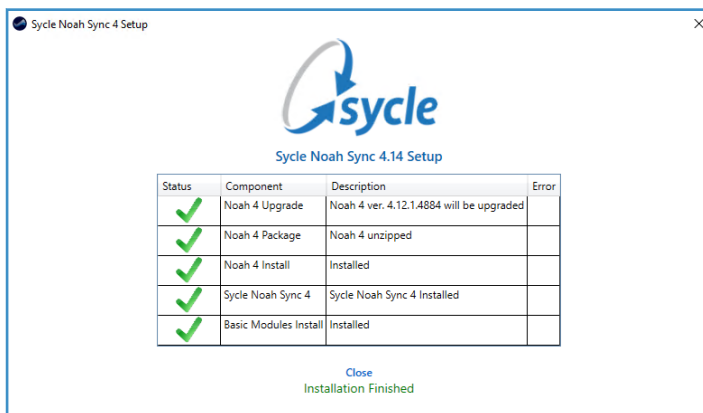

15. (If applicable) Un-pause or enable any non-Microsoft antivirus programs or firewalls that you disabled in [Step 1](#).
16. (If applicable) Repeat this procedure for all other PCs running SNS.

Appendix A. Glossary

The following terms are used in this document:

Sycle Customer Experience Team	<p>Sycle's dedicated team is just a phone call or email away.</p> <p>For sales, training or support requests, call 888.881.7925.</p> <p>You can also contact the team via email at <support@sycle.net>.</p>
Sycle Noah Sync (SNS)	<p>Sycle's SNS offers the ability to store your Noah data securely in the cloud, letting you share data between workstations and locations with no special networking configuration required. A patient's entire audiological record is available for view from any clinic that has access to the patient record.</p>

Appendix B. Document Revision History

Rev #	Date	Description
1.0	2023-04-17	<ul style="list-style-type: none"> • Changed directory structure in Creating a Noah System Data Backup, page 4 to include the entire /Noah (or Noah 4) directory. • Added Backing Up the Noah Directory, page 20 chapter to Upgrade section. • Added instruction in Sycle Noah Sync (SNS) Upgrade Requirements, page 19 to contact CE if current version lower than 4.9. • Added caution admonition to Running the Upgrade, page 22 for users to contact CE if current version lower than 4.9. • Added the upgrade backup instructions to the summary in Overview, page 1. • Document classification changed to Public. Removed Confidential label, <i>Legal Notice</i>, and watermark.
0.3	2023-04-12	<ul style="list-style-type: none"> • Minor changes to sentence structure in Sycle Noah Sync (SNS) Upgrade Requirements, page 19.
0.2	2023-04-04	<ul style="list-style-type: none"> • Added Overview, page 1 and the summaries of instructions. • Updated installation and upgrade requirements. • Changed download link for SNS installer package from sycle.com/pricing-features/noah-sync/ to sycle.com/noah-synch-download. • Added installation/upgrade steps for disabling and re-enabling any non-Microsoft anti-virus or firewall programs.
0.1	2023-03-31	Initial document creation.